

Air-polishing

ALL-IN-ONE FOR EASY PROCEDURES

Enhancing your workflows
through innovation

air-n-go
easy





Simplicity & efficacy all-in-one

For a supra- and a sub-gingival
treatment

A single handpiece

A single handpiece to cover a range of care
and treatment

AIR-N-GO® easy combines essential functions
with reliability to facilitate your everyday life:

- Fast start-up: direct connection to the high
speed handpiece connector
- Greater freedom of movement: 360° rotation,
fast and accurate treatment of areas
requiring care
- Easier cleaning of the handpiece: main parts
can be dismantled

Dual function

- Supra-gingival polishing **SUPRA
FUNCTION**: daily prophylactic polishing
- Sub-gingival polishing **PERIO FUNCTION**:
maintenance, periodontal and peri-implant
treatments

This medical device allies simplicity with
performance.

A complete solution

4 nozzles for active control of the
progress of periodontal disease



A full range of powders

- Sodium bicarbonate-based **CLASSIC** powders•
- Glycine-based **PERIO** powder


Classic & Perio



4 allies to prevent and fight periodontal disease

AIR-N-GO® nozzles, which result from the most advanced fluid dynamics modelling techniques, provide continuous, accurate and controlled spray whatever the clinical context.

In daily practice these nozzles provide power and efficiency for the removal of supra- and sub-gingival biofilm.

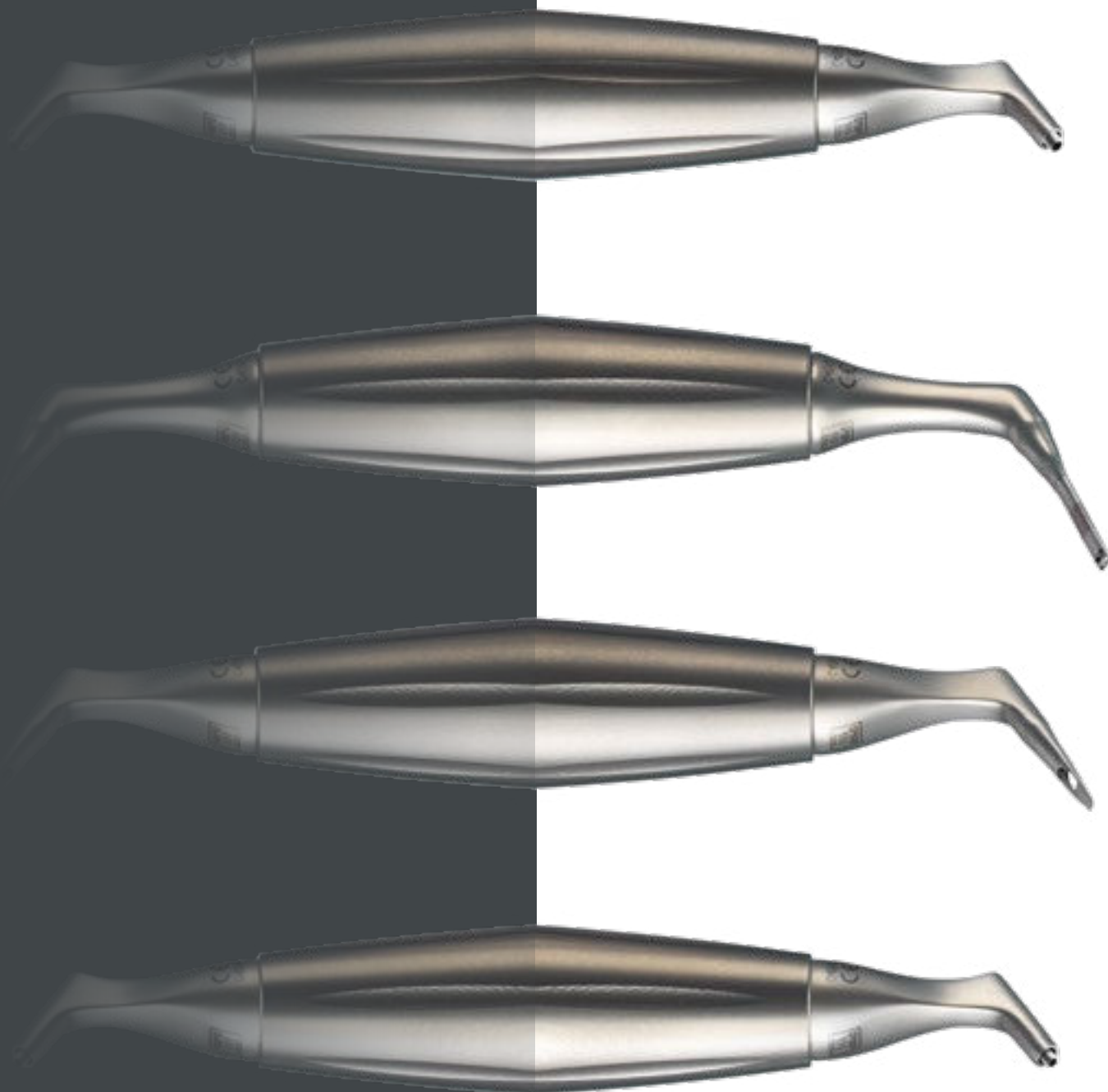

Air-N-Go® SUPRA nozzle



F10125 Supplied as standard

- Supra-gingival prophylactic treatment in combination with AIR-N-GO® CLASSIC powders (finishing after scaling, preparation for whitening, cleaning of inter-proximal areas, fissures and crevices, preparation of the tooth surface before etching, cleaning of orthodontic brackets, removal of supra-gingival biofilm, plaque, stains and temporary cement residues)
- Guarantees a fast, efficient procedure for greater working comfort



Dr P. ROUAS (France) clinical case

Air-N-Go® PERIO EASY nozzle



F10126 Supplied as an option

- Removal of the sub-gingival biofilm
- A real asset for non-surgical, sub-gingival treatment of teeth and implants (3-8mm periodontal pockets)
- Innovative design in a "shoehorn" shape providing easy sub-gingival insertion
- 5 seconds per face is all that is needed



Dr C. LIM (U.K.) clinical case

Air-N-Go® PERIO nozzle



F10127 Supplied as an option

- Sub-gingival use for surgical and non-surgical treatments (tunnelling and open flaps)
- Decontamination of deep periodontal pockets (8-10mm)
- Powerful double lateral spray



Dr A. PATEL (U.K.) clinical case

Air-N-Go® PERIO MAINTENANCE nozzle



F10127 Supplied as an option

- Periodontal maintenance (Pockets up to 4mm)
- Treatment of mucositis and return to a healthy clinical situation
- Reduced consumption of AIR-N-GO® PERIO fine powder particles in the maintenance phase
- Action against the progression of periodontal disease



Dr C. LIM (U.K.) clinical case

Customized solutions for your clinical needs

AIR-N-GO® powders have been specially designed and studied to produce effective treatment that is gentle on all treated surfaces. All these powders are compatible with the AIR-N-GO® easy, no change of handpiece is required.



5 fresh flavours based on 100% natural aromas or essential oils:
Bottles (4 x 250g): F10254 Neutral, F10253 Cola, F10252 Raspberry, F10250 Peppermint and F10251 Lemon


- Complete range of active **supra-gingival powders** which meet any patients' preferences
- Action of removing coloured stains and bacterial plaque, even in areas that are difficult to access
- Scrupulous cleaning giving instant results
- Better patient comfort: gentle and efficient with limited bleeding during treatment

Particles with softened geometry (x1000)



Treated enamel: surface preserved without impacts or scratches (x300)

Stratified particles with aggressive geometry (x1000)

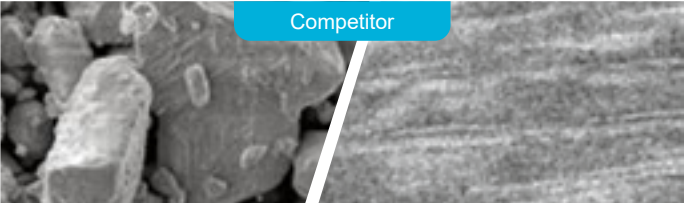

Treated enamel shows numerous scratches (x300)

Source: ICMCB CNRS (Institute of Condensed Matter Chemistry - French National Centre for Scientific Research) S.E.M. photos


air-n-go perio

Natural glycine-based (≈ 25 µm)

- Exclusive use for **sub-gingival treatments** as a complement to the initial therapy
- Low abrasiveness preserving the treated surfaces (teeth and implants) and the most fragile anatomical areas
- Prevents the progression of the sub-gingival biofilm
- **Multiple advantages with unique properties:**
 - Plaque removal: during periodontal maintenance on exposed root surfaces without sacrificing substance
 - Reduction of the bacterial load in the periodontal pockets



F10270 (3 x 160 g)

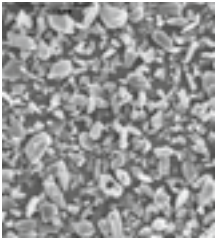

AIR-N-GO "PERIO" (x1000)

Source: ICMCB CNRS (Institute of Condensed Matter Chemistry - French National Centre for Scientific Research) S.E.M. photos



About Prophylaxis workflow

This combination of solutions helps dentists easily prevent and diagnose potential dental problems such as periodontal disease and tooth decay. It also allows to perform safe and minimally invasive procedures, providing best care and comfort for patients.



Diagnosis

To ensure a reliable diagnosis, I need to get a quick and clear view of initial signs of caries, periodontal disease, and other conditions.

C50



C50 intraoral camera is the perfect companion for non-invasive and reliable caries detection, thanks to its fluorescent technology.



Treatment

I wish to carry out fast and accurate treatment, while reducing hand fatigue.

newtron



bliss

Newtron® P5XS B.LED and Bliss™ is the perfect combo to preserve tissues with efficacy and precision. Less force is needed for more comfort.



Air polishing

It is important for me to control the progression of the disease, without neglecting the comfort of my patient.

**air-n-go
easy**



With Air-N-Go® easy, provide a gentle and effective treatment that limits bleeding. Stains and bacterial plaque are removed even in hard-to-reach areas.



Follow-up

In order to educate my patients for better treatment acceptance, I need to highlight the results of my dental care.

C50



C50 provides a reliable and immediate check of the clinical condition. The images can be shown to the patient and saved to his file.



Protection

I need to ensure optimal clinical protection for my patients and limit the risk of cross-contamination.



riskontrol

Riskontrol® disposable tips are your must-have solution to perform dental care in a safe and clean place. It assures protection, comfort and safety in every daily clinical procedures.



Keep optimal performances in a long-term



Easy monitoring and routine maintenance

Preventive cleaning after each use allows to maintain the Air-N-Go® easy performances.

Air-N-Go® Classic powder is cleaned with 4% acetic acid aqueous solution such as $C_2H_4O_2$ molecular formula spirit vinegar or diluted lemon juice.

Air-N-Go® Perio powder is cleaned with C_2H_5O formula 27% ethanol aqueous solution such as green Listérine®.


Easy disinfection

Clean the outside of Air-N-Go® easy with an alcohol disinfection wipe.



Easy sterilization

Air-N-Go® easy body and its nozzles are autoclavable.

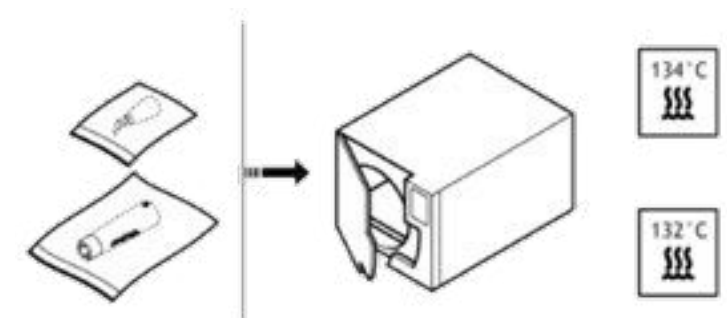



SATELEC

A company of ACTEON® Group

**17 avenue Gustave Eiffel,
ZI du Phare,
33700 Mérignac, France**

Tel. +33 (0) 556 340 607

Fax. +33 (0) 556 349 292

info@acteongroup.com

Follow-us

