MCDental.

Dental Decontamination.

ISSUE: 3





To ensure a safe and hygienic environment:

Covering a comprehensive range of supplies and accessories to ensure effective sterilization and disinfection.

AUTOCLAVES





Lyla S type Sterilizer

Your customised solution – now and in the future.

With its innovative duo-chamber system, lubrication in record time and the environmentally friendly Care Set, it impresses all along the line.

Overall dimension: 465 × 452 × 646 mm

Weight: 40 kg, 42 kg

Working range: 7 to 10 cycles

Lyla 17/22 loads: maximum solid 4 kg/typical 2

kg/maximum porous 0.5 kg

Lyla 17/22 container load: maximum 4 kg

SKU: 19951154, 19952154





Lara B type Sterilizer

Incredible Inside & Out - As unique as your practice.

Built for highest demands: incredibly user-friendly thanks to EliSense information system via LED indicators and display.

Overall dimension: 465 x 452 x 646 mm

Weight: 42.5 kg, 44 kg

Working range: 7 to 15 cycles

Lara 17 loads: maximum solid 4.5 kg/typical 2

kg/maximum porous 1.5 kg

Lisa 22 loads: maximum solid 5.5 kg/typical 2

kg/maximum porous 2 kg

£6,857 £4,65

£7,700 £5,2

SKU: 19941120, 19942120, 19933110





Price Includes:

Lisa Mini Sterilizer

Speeds up your daily workflow.

Reprocessing of unwrapped loads in less than 10 minutes, and much more. The fully equipped type B sterilizer offers the fastest cycles in its performance class.

- Simple: the home screen immediately displays the cycle selection when switched on.
- · Risk-Free: thanks to the possibility of customizing the start screen with preferred cycles, daily activities are simplified.

5 L Sterilsier features integrated vacuum pump

£4,900

SKU: 19945100





Lisa B Type Sterilizer

Incredible Inside & Out - As unique as your practice.

The Lisa sterilizers – built for highest demands: incredibly user-friendly thanks to EliSense information system via LED indicators and display.

Overall dimension: 465 x 452 x 646 mm

Weight: 46 kg, 47.5 kg

Working range: 7 to 15 cycles

Lara 17 loads: maximum 4.5 kg/typical 2

kg/maximum porous 1.5 kg

Lisa 22 loads: maximum solid 6 kg/typical 2

kg/maximum porous 2 kg

£8,467 £5,5

£9,580 £6,2

SKU: 19931184, 19932184

Dental Decontamination Supplies: LUBRICATION SYSTEMS





Assisting Twin

Increase your handpieces lifetime and reduce the costs of repairs!

With its innovative duo-chamber system, lubrication in record time and the environmentally friendly Care Set, it impresses all along the line.

W&H Assistina Twin is ergonomically designed to work effectively and efficiently making it even quicker and easier to maintain your handpieces through the innovative oil nebulization technology.

£3,202.00

SKU: 30310003





Assistina One

The ONE that speeds up MY daily operation

Fast, thorough and easy to use: Assistina One speeds up daily operations and facilitates dental handpiece care. With the shortest cycle time in its class and its high-quality care results, the innovative maintenance solution is a top dental device designed for top dental professionals.

Entire lubrication and spray channel cleaning process only takes 18 seconds

£2,169.00

SKU: 30398003



quick

CONNECT by W&H



Assisting Twin Care Set

W&H Assistina Twin Care Set to ensure the correct functioning of the Assistina Twin.



- · Maintenance for your instruments safely and quickly.
- Its productivity is equivalent to 20 bottles of spray oil.
- With a fill level of 400 ml per spray oil bottle.

SKU: 08064970

Assisting One Care Set

Handpiece care with the W&H Assistina One Care Set is more cost-effective than manual lubrication or other systems.

Product Includes:

- 1 cartridge of Activefluid (200 mI)
- · 1 cartridge of Service Oil F1 (200 ml)
- 1 HEPA filter



SKU: 08064970






Autoclave Dust Filter

- · Large dust filter for Lisa pre-2019.
- Lina & MS-22 all models
- 400 cycles to protect internal components



SKU: F364502X

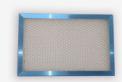

Dust Filter Small

- Lisa/Lara*/Lyla current models.
- (LARA XL does not have a dust filter)
- 400 Cycle frequency.

SKU: F364511X



Door Seal

• W&H Door Seal for 800 Cycle frequency for Lisa/Lara/Lyla models.



SKU: F460504X



Bacterial Filter

For all W&H Sterilizers.

- Easy to install, no tools required.
- To prevent the penetration of the bacteria in the chamber during the sterilization process.

SKU: W322400X

PRICES AND PRODUCT AVAILABILITY ARE SUBJECT TO CHANGE WITHOUT PRIOR NOTICE. VAT NOT SHOWN.

ULTRASONIC CLEANERS





MHC Ultrasonic Cleaner

- AC 230V, 50Hz / AC 110V, 60Hz.
- Frequency 40 KHz.
- Transformation of the ultrasounds in mechanical vibrations by piezoelectric transductors.
- Programmable from 1 at 99 minutes by step of 1 minute.
- Display of the remaining time.

Without heating





SKU: MHC30, MHC50





MHC Ultrasonic Cleaner

- AC 230V, 50Hz / AC 110V, 60Hz.
- Frequency 40 KHz.
- Transformation of the ultrasounds in mechanical vibrations by piezoelectric transductors.
- Programmable from 1 at 99 minutes by step of 1
- · Display of the remaining time.
- Heating activation / interruption regardless of the cleaning cycle
- Programmable from 20 to 80 ° C

With heating



SKU: MHC30C, MHC50C



DDL4 Super Jug Water Distiller

- Produces 4 Litres of de-ionized water on demand in 20 minutes!
- No electricity supply or plumbing needed & the unit is portable.
- Includes tower water dispenser driven by a cell battery for dispensing metered de-ionized water to autoclaves and DUWL bottles.
- Includes TWO additional cartridges in the delivered package offering de ionized water for up to 18 months.



SKU: DDL4





MHC Water Distiller

The MHC water distiller MHC produces perfectly pure water by distillation. The machine heats the water up to 100° which rises up as steam to the condenser.

Power supply: AC 230V - 50Hz

Power: 750W Tank capacity: 4L

Speed distillation: 0,7 to 1 L/h Safety Temperature: 160°

Weight: 4 Kg



SKU: MHC400

General Lubrication Sys

Maintenance for after use and before











- Detach handpiece from tubing.
- Remove bur from Handpiece.
- DO NOT submerge handpiece in any solution.



STEP: 2 CLEAN

- Wipe away any external debris using a non-alcoholic approved disinfectant wipe.
- Using a dry tissue, remove any excess detergent leftover on the handpeice from the wipe.
 - DO NOT use any chemical solution or submerge in water.



STEP: 3 CONNECT

- Select the appropriate motor for the handpiece model you are lubricating.
- Connect the motor securely to the lubrication system, ensuring a tight fit to prevent leaks during lubrication phase.

Make sure to point the handpiece heads towards the absorbent pads whenever available. If absorbent pads are not included as a standard part of the system, using clean paper towels for each cycle would suffice.

Best practices for cleaning and maintaining Equipment:

Avoid Chemical Disinfectants: Refrain from using chemical disinfectants, abrasive cleaners, or antibacterial

Avoid Submerging: Avoid submerging any handpieces into any type of liquid for cleaning or disinfection.

Avoid Ultrasonic Cleaners: Do not utilize ultrasonic cleaners for cleaning any handpieces.

Ensure Proper Installation: Ensure all instruments are removed from your handpiece or equipment prior to c

Autoclave: When autoclaving, do not exceed a temperature of 138°C and do not overfill trays as this will prev

stem Guide.

Autoclaving.





STEP: 4 POWER

- Power on the lubrication system according to the manufacturer's instructions.
- Allow the lubrication system to run its cycle, ensuring that it completes the lubrication process thoroughly.
- Wait until the lubrication cycle finishes before proceeding to the next step.
- Check absorbent pads or paper towels for any discoloration indicating oil contamination. If discolorations observed, promptly seek guidance from MC Repairs.





STEP: 5

- After lubricating the handpiece, remove it from the adaptor.
- Seal the pouch according to the manufacturer's instructions, ensuring a tight seal to maintain sterility.





STEP: 6 STERILISE

- Before placing your equipment in the autoclave, ensure you've selected the correct cycle and verified that the autoclave temperature does not exceed 138°C
- Load the items to be sterilized into the autoclave chamber according to the manufacturer's guidelines.
- For any further advice or guidance, don't hesitate to contact

Please Note: these are generalised guidelines as recommended by MC Repairs Ltd. For exact guidelines, please refer to Lubrication System user manual.

soaps as cleaners.

MCRepairs.

leaning and lubricating.

ent full sterilization.

AUTOCLAVES





NSK IMax S S-Type Autoclave

The NSK iMax S S-Type Autoclave provides all the benefits of a traditional displacement autoclave but is effective for the correct sterilisation of handpieces.

20 min fast cycles for unwrapped loads. 25 min handpiece cycle Affordable, effective upgrade from non-vacuum autoclaves Full vacuum drying phase (Dry handpeices!)

External dimensions: 60cm x 45cm x 43cm Max Load: 4kg (solid), 1.5kg (porous)



SKU: DX91803S0200





NSK IClave Plus Autoclave

The panel's sophisticated design and excellent visibility make it easy to operate and maintain, with color changes enhancing sterilization reliability.

Advanced heating system leveraging excellent thermal conductivity. NSK's innovative heating system optimises the high heat conductivity of copper.

External dimensions: 445 x 532 x 428 mm Maximum Power Consumption: 1,900W Air Expulsion System: Vacuum pump 1,3,4 vacuum Max Load: 4kg (solid), 1.5kg (porous)



SKU: DXY1003077





iClave 24 Autoclave

The new iClave 24 autoclave perfectly combines performance, user experience, and technology. Ensuring a comfortable and worry-free experience.

In addition, the new high-efficiency steam generator allows for superior performance, creating perfectly sterilised equipment in less time.

External Dimensions: W 451 x H 492 x D 669 mm Cycle Water Consumption*: 700 ~ 800 cc Tank Capacity (clean and used water): 5 litres

Sound Level: 60 dB(A)



SKU: 92412B0100



NSK Printer





Helix Test

Allows you to verify steam penetration into hollow instruments.

This test kit consists of a plastic tube with an indicator strip housing on one end and is supplied with 250 easy-to-read test strips compliant with the EN 867-5 standard.



SKU: 9900009

SKU: 9900059

Dental Decontamination Supplies: LUBRICATION SYSTEMS





ICare Handpiece Oiling Unit

NSK iCare Handpiece Oiling Unit sets a new standard for handpiece maintenance systems. In one stylish and compact body, it offers you both ease of use and multiple functions to achieve complete cleaning and lubrication.

Power Supply: AC 100-240V 50/60Hz Air Pressure: 0.5-0.6 MPa (5.0-6.0 kgf/cm2) Dimension: W280 x D240 x H350 (mm) Solution Tank Capacity: 1.2L





SKU: Y1002787, Y1002788





ICare+ Decontamination Unit

Manual cleaning of instruments does not remove all deposits and residues. Dynamic instruments that are poorly maintained or not adequately cleaned may become damaged and pass on infections to both patients and practitioners.

Any instrument used in dentistry must be perfectly cleaned, disinfected and maintained. NSK is aware of the requirements in the daily practice of dentistry and, after many years of research, it has come up with a simple and effective solution thanks to iCare:

£3,995

SKU: S102001, S103001



PTL Adaptor

B2/M4 Adaptor

KV Adaptor

SR Adaptor

WH Adaptor

BA Adaptor

F-Type Adaptor



iCare/iCare+ Adaptors

• NSK iCare Handpiece Adaptors for use with the NSK icare and icare C3 machine.

FROM

£110.00

SKU: Z1127010, Z1127011, Z1127013 Z1127012, Z1127014

Mist Filter Set

• Pack of 12, for iCare/Care3



SKU: Y900159



HELIX Test

- · Helix Test allows you to verify steam penetration into hollow instruments.
- Includes a tester and 250 indicators

£184

SKU: DX9900059



- · Removes limescale and dirt resi-
- due from the chamber.
- Box of 10

SKU: DX0230050





iCare+ nClean

- Cleaning solution for NSK iCare Plus Maintenance System.
- 6 x 500ml Bottles



SKU: ACL600



Sealer Newseal

· With high output, easily and automatically seals envelopes before sterilization.



SKU: DX100071



iCare+ nCid

- · Disinfection solution for use in NSK iCare+ Maintenane Unit.
- 6 x 500ml Bottles.



SKU: ACD600



Purity for iClave plus

· Fits on top of the iCare plus automatically filtering water from the mains supply.

SKU: DX98018030

AUTOCLAVES





Mocom B Futura Autoclave

Designed to make your life easier, reduce costs and guarantee higher efficiency. Equipped with advanced connectivity and traceability systems.

The system boosts surgery efficiency by resolving demineralized water management problems, eliminating storage issues and component deterioration due to poor water quality.

External dimensions: 480 x 500 x 600 mm 3 year or 3,000 cycle Warranty 12 month tap water filter

£4.950

£5,450

£5,950

SKU: bfutura17, bfutura22, bfutura28





Mocom B Classic Autoclave

The Mocom B Classic Autoclave makes work easier by optimising times and costs.

High-level quality and functions guarantee safety, reliability and ease of use. Three versions available: 17, 22 and 28 litres. Thanks to the same external dimensions and LCD display, B Classic is a perfect solution for everyone.

3 year or 3,000 cycle warranty

External dimensions: 445 x 532 x 428 mm

£3,950

£4,450

£4,950

SKU: bclassic17, bclassic22, bclassic28





Mocom Supreme Autoclave

A closed-loop sterilizer with an internal filtration system recycles water for reuse, reducing waste, costs, and increasing efficiency in dentistry surgeries. It resolves demineralized water storage and supply issues.

50 sterilization cycles with just 5 litres of tap water.

External dimensions: 480 x 500 x 600 mm 3 year or 3,000 cycle Warranty 12 month tap water filter

£5.950

£6,450

£6,950

SKU: supreme17, supreme22 supreme28





Mocom Tethys H10 Plus

The Evolution of Disinfection

Designed to give you the best. The best performance in a simple and fast work-flow. An innovative device that replaces the numerous manual tasks typical of the stages preceding sterilization, thus reducing the operator's workloads. Tethys H10 Plus makes the instrument reconditioning process even simpler and more functional.

3 year or 3,000 cycle warranty

Price includes: Installation, Training, Comissioning.

£4,500

SKU: TethysH10

LUBRICATION SYSTEMS







Kavo QUATTROcare PLUS 2124A

The cost-efficient solution for optimum care and long-term maintenance of value.

KaVo QUATTROcare maintenance system, KaVo sets standards of efficiency for handpieces care. The complete care procedure for four handpieces, only takes one minute

- Very easy to use
- Four handpieces in one minute
- High efficiency
- Integrated chuck care program

£2.501.00

SKU: 1.008.3805

6 FREE cans with every **QUATTROCARE order!**





QUATTROcare Plus Spray

- · Designed to clean, lubricate and purge your instruments with the simple push of a button.
- KaVo's QUATTROcare Plus Spray cleaning and care of instruments and air motors is easy. KaVo's QUAT-TROcare spray is designed specifically to withstand sterilisation.
- •1x500ml can
- Environmentally friendly propellent

1 x 500ml £40.00

6 x 500ml

SKU: 1.005.4525

QUATTROcare Adaptors

KaVo Quattrocare Adaptors & Chuck for 2104, 2104 A, 2124 A.

 KaVo QUATROcare Adaptors connects the handpieces to the QUATTROcare maintenance system to be cleaned, purged and lubricated.



SKU: 0.411.7603

